

Operative



PRECLINICAL GUIDE



Amalgam

Tech Operative

Amalgam : An alloy of mercury with one or more other metals.



Dental Amalgam Alloy

- Contains silver (Ag), tin (Sn), copper (Cu), zinc (Zn), palladium (Pd), indium (In)

Dental Amalgam

- Formed by mixing alloy particles with mercury in the dental office to form a **plastic mass**
- The mass is inserted into the prepared cavity while soft, then allowed to **harden intraorally**

Properties of Amalgam

Advantages .

Mechanical Properties

- High **compressive strength**
- High **fracture resistance** under occlusal forces
- Adequate **form stability** (maintains anatomical form and proximal contacts)
- Insoluble, **high wear resistance**, and **low creep value** in recent types

Biological Effects

- Self-sealing** over time due to corrosion products
- Microleakage, sensitivity, and recurrent caries may occur initially

Other Properties

- Low cost**, low technique sensitivity, easy manipulation
- Widely used with **high durability**

Dis- Advantages .

- Mechanical** : Low **tensile strength** (~25% of compressive strength) → brittle

Requires **sufficient thickness**; prone to fracture if thin (marginal ditching)

- Physical** : High **thermal conductivity** → requires protective base in deep cavities
- Aesthetically** : Objectionable **metallic color**, tarnish (**discoloration**) + **corrosion** (**destructive breakdown of the metal**)





Tarnish and Corrosion


Mechanical Flow & Creep

- **Flow:** deformation before setting
- **Creep:** time-dependent plastic deformation after setting under constant load


Consequences:

- Marginal ditching
- Flattening of contact
- Gingival overhangs
- Saucering of occlusal anatomy

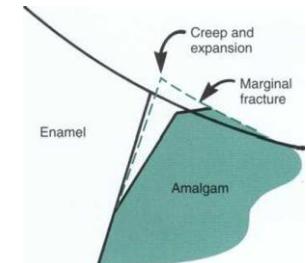

Elastic Modulus

- High modulus → rigid restoration
- Causes stress concentration and poor stress distribution


Galvanic Corrosion

- Occurs between dissimilar metals in presence of moisture 
- Leads to current flow and corrosion 
- Causes **post-operative sensitivity** 


Uses of Amalgam.


- Posterior teeth due to metallic color 
- Small Class I and II restorations with sufficient sound tooth structure 
- Core build-up under full coverage restorations 

 **Note:** Composite resin increasingly replaces amalgam due to material improvements



Contraindications

1. Extensive lesions with undermined cusps (adhesive restorations preferred)
2. Esthetic zones 
3. Presence of opposing metal restorations (to avoid galvanic activity) 



Cavity Preparation for Amalgam

Based on the brittle and non-adhesive nature of amalgam, cavity preparation must:

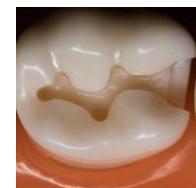

1. Be **conservative** 
2. Provide **adequate bulk** 
3. Place amalgam in **self-cleansable areas** 
4. Ensure **supported tooth structure** 

Cavity Outline Form

- **Depth:** 1.5–2 mm or 0.5–1 mm beyond DEJ 
- **Width:**
 - Buccolingual: 1/4–1/3 intercuspal distance
 - Mesiodistal: include triangular fossa (minimum marginal ridge thickness 1.5 mm)
- **Wall direction:**
 - Buccal & lingual walls slightly convergent
 - Mesial & distal walls slightly divergent
- **Outline:** Sweeping curves 
- **Finishing walls:** Smooth and rounded; flat, uniform floor 

Class II:

- Proximal flaring to open contact
- Isthmus wider than cavity to avoid necking of amalgam
- Proximal walls follow external tooth contour
- CSA= 9

Composition of Amalgam Alloy

Element	Strength	Setting Expansion	Setting Time	Creep	Tarnish & Corrosion Resistance
Ag					
Sn					
Cu			—		
Hg					
Zn	<ul style="list-style-type: none"> • Acts as scavenger to prevent oxides formation during alloy manufacturing • Improve plasticity and workability of amalgam 				
In, Pd	Enhance mechanical properties				





Classification of Amalgam Alloys

1 According to Zinc (Zn) Content

● Zinc-containing alloys = $Zn > 0.1\%$

○ Zinc-free alloys = $Zn < 0.1\%$

- Zinc acts as a scavenger during alloy manufacturing ↗
- Prevents oxide contamination during manufacturing
- ⚠ If moisture contamination occurs:

→ Excessive delayed expansion → Due to formation of hydrogen gas (H_2)

the ↑↑ Cu amalgam the most corrodible Y₂ phase is eliminated
But if excess Hg is present → 2γ phase

2 According to Copper Content


Type	Copper %	Composition	Particle Type	γ ₂ Phase
Conventional (Low Cu alloy)	4% Cu	Ag–Sn alloy	Lathe-cut	Present
High Copper Alloys				
Admixed High Cu Alloy	9–20%	Mixture of: • Low Cu alloy particles • Ag–Cu eutectic alloy • Eutectic alloy: 72% Ag – 28% Cu	Lathe-cut + spherical	Eliminated
Unicompositional High Cu Alloy	13–30%	Single type of alloy	Spherical	Eliminated

3 According to Particle Shape


- Lathe-cut alloys
- Spherical alloys ○
- Admixed alloys

Based on shape of alloy particles


Lathe cut Alloys
Irregular shaped particles

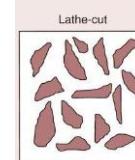

Spherical Alloys
Regularly shaped particles

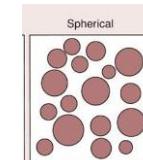

Admixed Alloys
Lathe cut + Spherical

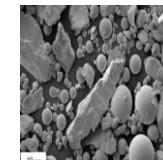

4 According to Particle Size

- Coarse-cut
- Fine-cut
- Micro-cut



Properties of Mercury



- Only metal liquid at room temperature ↗
- High surface tension
- Combines readily with Au, Ag, Sn, Cu, Zn
- **Highly toxic** ☠






Mercury Toxicity



- Dental staff are at greater risk than patients 
- Toxic manifestations: paresthesia, ataxia, joint pain, hearing loss, death



Main Sources of Exposure

- Accidental spills
- Poor mercury hygiene (exposure, ventilation, storage)
- Direct contact
- Amalgamator
- Removal of old restorations

Most serious hazard: Mercury vapor 





Mercury Hygiene Procedures

- Use non-touch techniques  
- Mechanical amalgamators with good seals
- Store mercury under sodium thiosulfate solution 
- Immediate spill cleanup and proper ventilation 
- Use pre-capsulated amalgam 




Reaction of Amalgam Alloy with Mercury



Low Copper Alloys + Mercury



Reaction



- Mercury **should not be enough** to react with all γ phase




Reaction Stages


- In the beginning of formation of γ_1 & γ_2
- the mix is **plastic** → can be **condensed** into the cavity
- With time → γ_1 & γ_2 crystals grow → amalgam **hardens** 

Conventional Low Copper amalgam


Set material: Unreacted particles (γ) surrounded by a matrix of the reaction products (γ_1 & γ_2).

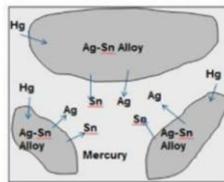

Gamma (γ) = Ag_3Sn
unreacted alloy
strongest phase and corrodes the least
forms 30% of volume of set amalgam

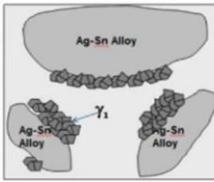

Gamma 1 (γ_1) = Ag_2Hg_3
matrix for unreacted alloy and **2nd strongest phase**
10 micron grains binding gamma (γ) 60% of volume

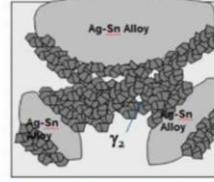


Gamma 2 (γ_2) = Sn_8Hg
weakest and softest phase
corrodes fast, voids form corrosion yields Hg which reacts with more gamma (γ) 10% of volume which decreases with time due to corrosion


High Copper amalgam

Set material: Unreacted particles (γ) and unreacted Ag-Cu surrounded by a halo of Cu_6Sn_5 (η) embedded in a matrix of (γ_1).

Admixed...Initial stage **Admixed...Final stage**

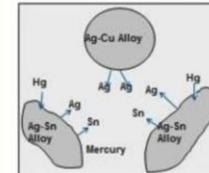

Ag enters Hg from Ag-Cu spherical eutectic particles.
Both Ag and Sn enter Hg from Ag_3Sn particles

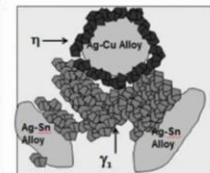


Sn diffuses to surface of Ag-Cu particles reacts with Cu to form (η) Cu_6Sn_5 (η) around unreacted Ag-Cu particles

High Copper Alloys + Mercury

1 Admixed High Cu Alloys + Hg

Alloys composed of TWO different particles:

- ◆ (A) Conventional Lathe-cut Low Cu Particles Same reaction as low Cu alloy:

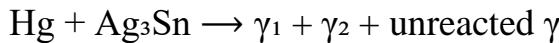

- ◆ (B) Ag-Cu Eutectic (Spherical Particles) Solid-state reaction begins:



🧠 Summary (Admixed)

- **γ_2 phase:**
 - Formed around **low Cu particles**
 - **Eliminated around eutectic particles**
- Net result:
→ Much less γ_2 than low Cu alloys

🧪 Solubility of Hg

Hg dissolves in metals in this order:

Sn > Ag > Cu

Metal	Hg solubility
Sn	170 mg
Ag	10 mg
Cu	1 mg

📌 That's why:

- Reaction of **low Cu alloy with Hg** happens first
- Then reaction with **Ag-Cu eutectic**





2 Unicompositional High Cu Alloys + Hg

Composition


- One type of particles containing:
Ag + Cu + Sn

Reaction Characteristics

- Elimination of γ_2 phase in ONE step
- Due to difference in Hg solubility:
Sn > Ag > Cu

Reaction Sequence

- Sn reacts first with Hg at particle periphery → temporary γ_2 formation
- Sn decreases → formation of **Ag–Cu eutectic**
- Final reaction:

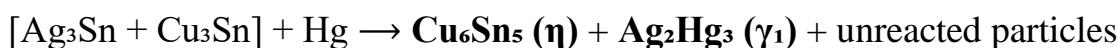

Summary (Unicompositional)

- Same particles:
 - Act first as γ phase
 - Then act as **eutectic**
- Final structure:
 - $\eta + \gamma_1$
 - **No γ_2**


Microstructure of Amalgam

1 Low Copper Amalgam

- Unreacted γ cores
- Held together by:
 - γ_1
 - γ_2 (matrix)
- Voids present
- Weakest structure

2 High Copper Amalgam

Feature	Admixed High Cu Alloys	Unicompositional High Cu Alloys
Alloy particles	γ phase & eutectic phases as unconsumed alloy particles	Unconsumed alloy particles
Unreacted cores	Present	Present
Surrounding structure	Unreacted cores surrounded by meshes & rods of η phase	Unreacted particles surrounded by meshes & rods of η phase
Binding phase	Held together by γ phase	Held together by γ_1 phase
Eutectic phase	Eutectic particles surrounded by η layers	
Matrix	γ_1 matrix, may contain η crystals	

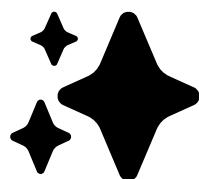

Thank You

